

*ENDOWRIST*<sup>®</sup>  
**STAPLER 30**  
FOR THE *DA VINCI XI*<sup>®</sup> SURGICAL SYSTEM



SURGEON CONTROL & AUTONOMY

FULLY WRISTED ARTICULATION

SMARTCLAMP FEEDBACK

## Ordering Information

### STAPLER 30 INSTRUMENTS

Stapler 30 shares the same cannula accessories and sheath as the *da Vinci Xi* Stapler 45.

**Straight Tip:** 470430



**Curved-Tip:** 470530



### STAPLER 30 RELOADS


All reloads include a new knife. All reloads have 6-rows of staples.



**Gray:** 48630M



**White:** 48630W

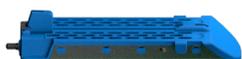

**Blue:** 48630B



**Green:** 48630G

# EndoWrist® Stapler 30 for the *da Vinci Xi*® Surgical System

The Stapler 30 Curved-Tip provides enhanced placement and visualization around tenuous vasculature. Combined with the *da Vinci Xi* Surgical System's crystal clear 3DHD vision, access, and surgeon control and autonomy, the *EndoWrist* Stapler 30 naturally expands the *da Vinci* stapling portfolio.

#### SURGEON CONTROL & AUTONOMY

Positioning, gripping, clamping, and firing of the *EndoWrist* Stapler are completely controlled through the *da Vinci* Surgeon Console.

#### FULLY WRISTED ARTICULATION

The *EndoWrist* Stapler has the greatest side-to-side articulation currently available, with a range of 108° total side-to-side. And unlike other staplers currently available, the *EndoWrist* Stapler articulates through a range of 54° total up-and-down, providing a fully-wristed cone of articulation.

#### SMARTCLAMP FEEDBACK

*SmartClamp* feedback reduces the guesswork associated with conventional stapling. Prior to firing, *SmartClamp* detects if the Stapler jaws are adequately closed on the tissue for the installed reload color.

#### Important Safety Information

Serious complications may occur in any surgery, including *da Vinci*® Surgery, up to and including death. Examples of serious or life-threatening complications, which may require prolonged and/or unexpected hospitalization and/or reoperation, include but are not limited to, one or more of the following: injury to tissues/organs, bleeding, infection and internal scarring that can cause long-lasting dysfunction/pain. Individual surgical results may vary. For Important Safety Information, indications for use, risks, full cautions and warnings, please also refer to [www.davincisurgery.com/safety](http://www.davincisurgery.com/safety) and [www.intuitivesurgical.com/safety](http://www.intuitivesurgical.com/safety).

The *EndoWrist*® Stapler 30 for the *da Vinci Xi*® System (IS4000) is not compatible for use with the *da Vinci*, *da Vinci S*, or *da Vinci Si* Surgical System. The *EndoWrist* Stapler 30, and *EndoWrist* Stapler 30 Reloads for the *da Vinci Xi* System are intended to be used with the *da Vinci Xi* Surgical System (IS4000) for resection, transection, and/or creation of anastomoses in General, Thoracic, Gynecologic and Urologic surgery. The device is indicated for adult and pediatric use. The device can be used with staple-line or tissue-buttressing materials.

The *EndoWrist* Stapler 30 and *EndoWrist* Stapler 30 Reloads for the *da Vinci Xi* System should not be used on tissue such as the liver or spleen, where tissue compressibility is such that clamping of the instrument would be destructive. Do not use the *EndoWrist* Stapler 30 or *EndoWrist* Stapler 30 Reloads for the *da Vinci Xi* System on the aorta.

© 2016 Intuitive Surgical, Inc. All rights reserved. Product names are trademarks or registered trademarks of their respective holders. PN1024875 Rev A 03/2016